



NUCLEO 64 CONNECTOR

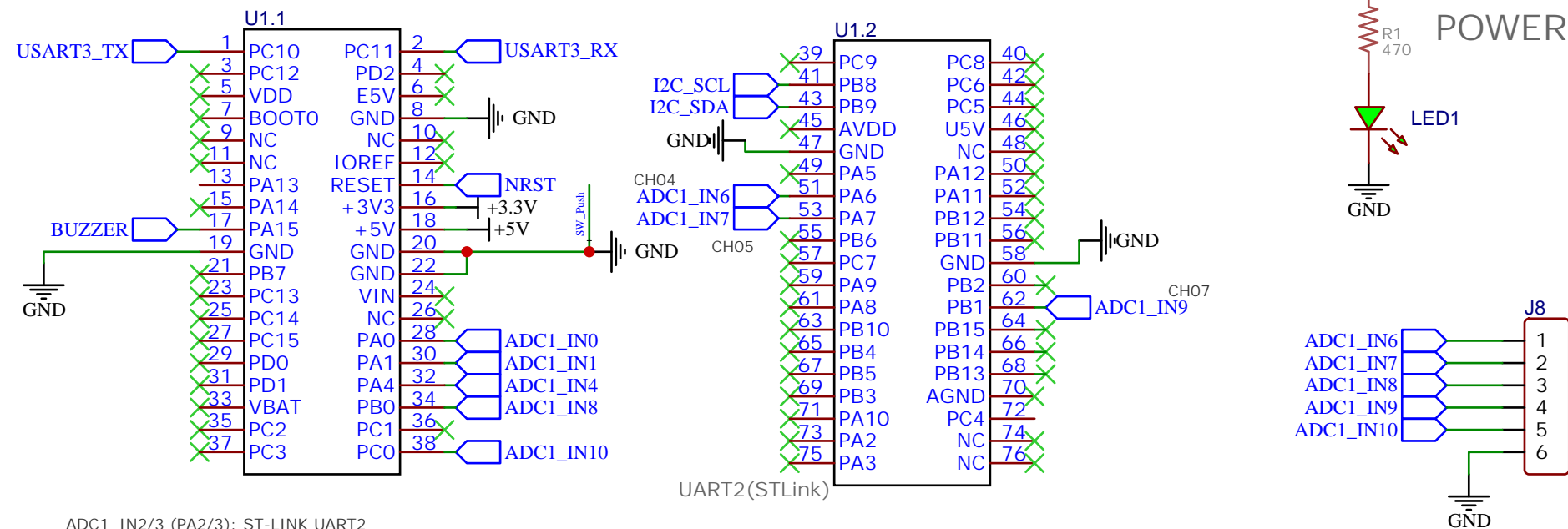


ADC1_IN2/3 (PA2/3): ST-LINK UART2
ADC1_IN5(PA5): LD2(Nucleo USR LED)

ADC1_IN0 : Ch01
ADC1_IN1 : Ch02
ADC1_IN4 : Ch03
ADC1_IN6 : Ch04
ADC1_IN7 : Ch05
ADC1_IN8 : Ch06
ADC1_IN9 : Ch07
ADC1_IN10 : Ch08



POWER

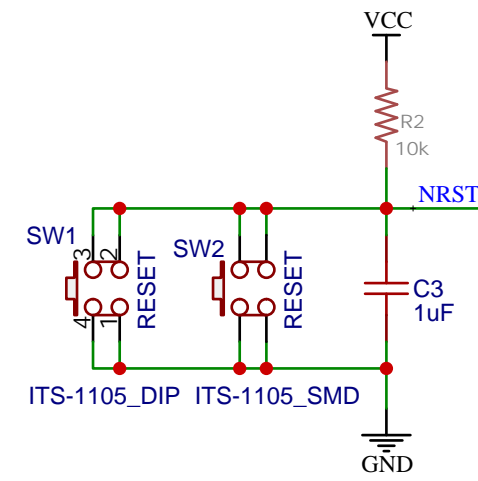

I2C

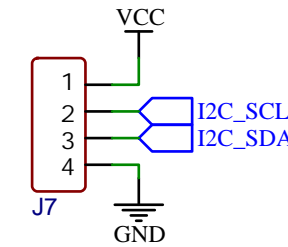


USART3

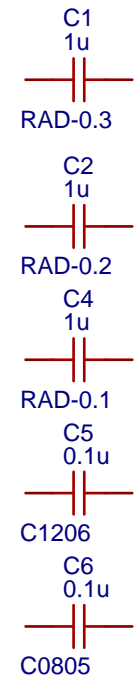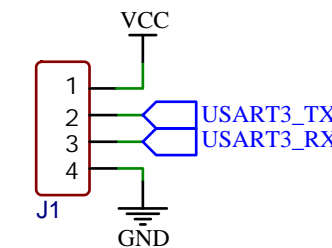


Sample Package Footprint

Power(Vcc)Selection

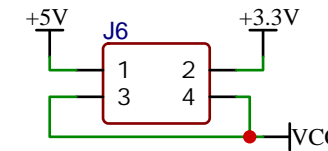


Electrode reference

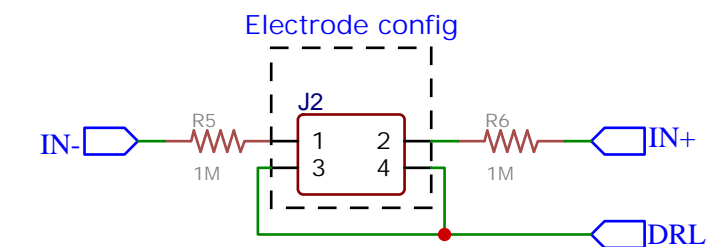


Bridge all three pads of the solder jumper to configure for 2 Electrode operations instead of the 3 electrodes!
Note: 2 electrodes will give noisier output than 3 electrodes!


Power Supply Noise Filtering

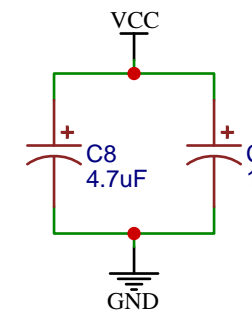

Amp Ref + Driven Right Leg (DRL)

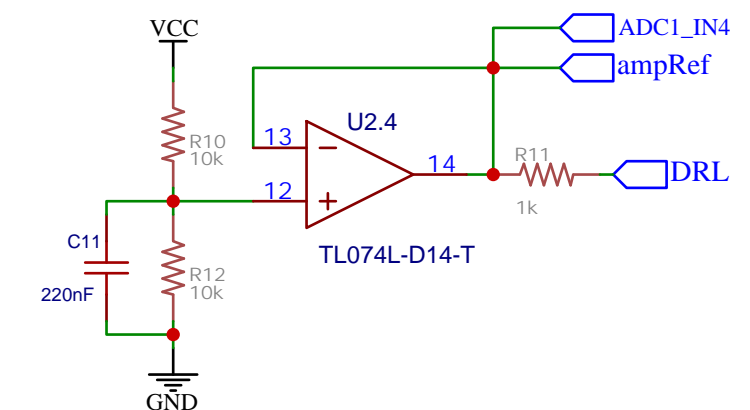


Instrumentation Amp (INA)

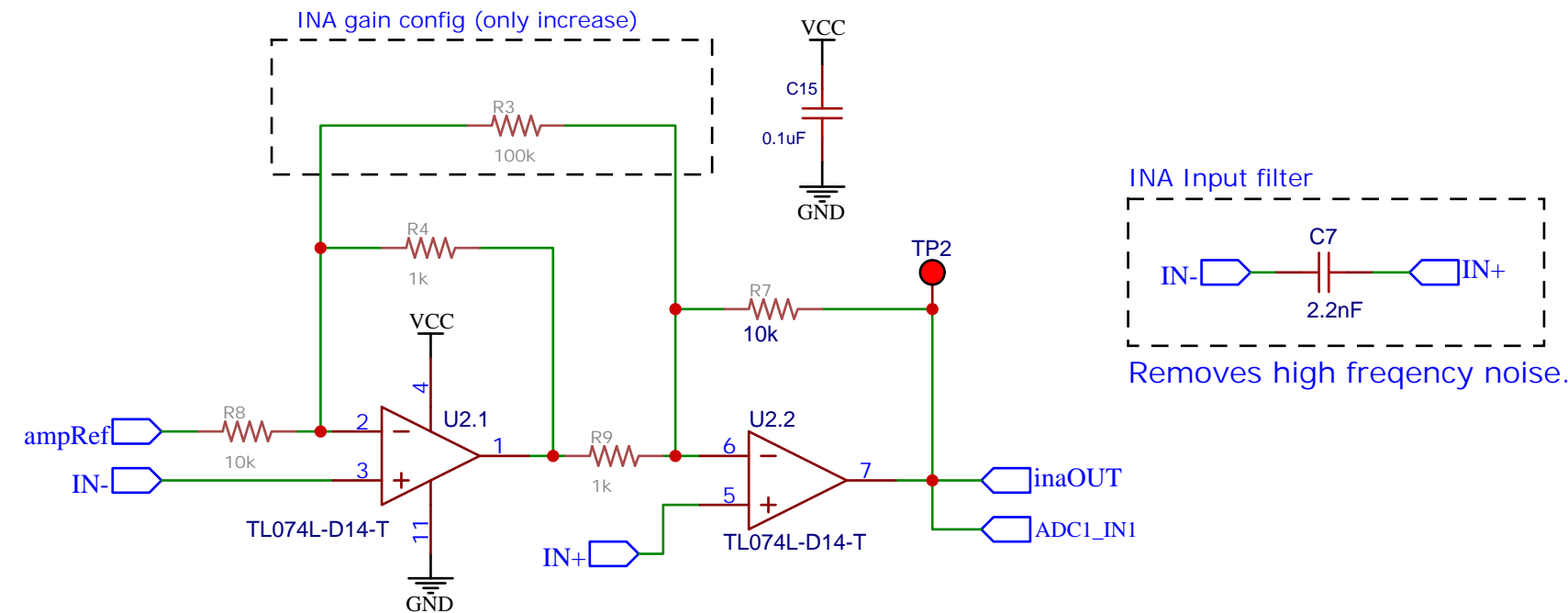


Create solder joint on J3 to narrow the Band Pass frequency range
Default: Wide input frequency band, use when recording EMG, EOG
Configured: Narrow input frequency band, use for EEG, EOG, and ECG

1000x Gain (Bandpass)

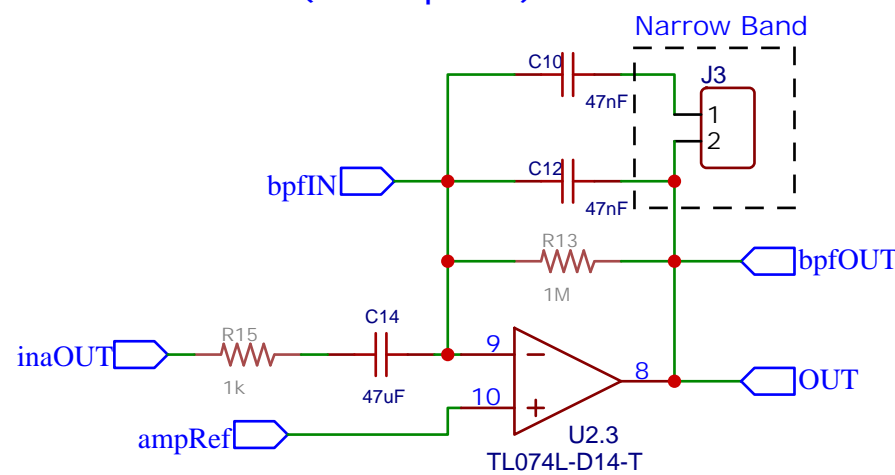


Band Pass Filter Configuration (BPFC)

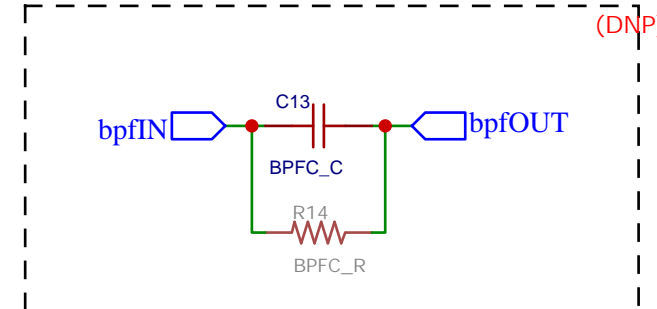

* C8 (BPFC_C) & R3 (BPFC_R) on the back side of the PCB can be used to configure the bandpass filter.
* Use them to configure output Gain (decrease) & Band (frequency range).

Header pins / Connectors

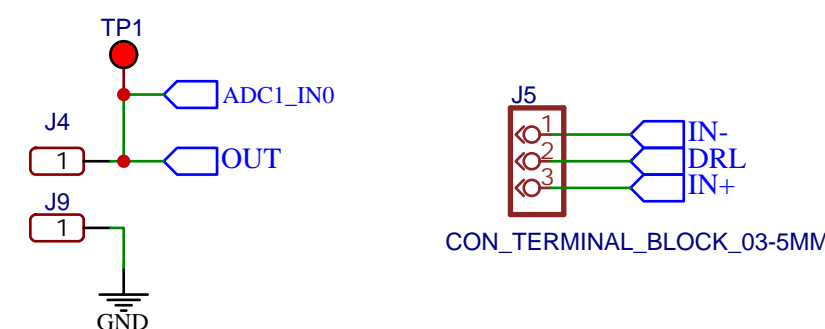

| | | |
|--|--|-------------------|
| TITLE: Bioamp ExG Experiment Brd.[KKU-BME] | | REV: 1.0 |
| Company: KKU-BIMSRL | | Sheet: 1/1 |
| Date: 2023-12-01 | | Drawn By: jwlee95 |